



A Hill-Rom Company

## HANDY LOW-FRICTION PRODUCTS



Always at Hand  
for Smooth Transfers




*Simple aids which are always at hand. This is how one could summarize Liko's new product series Handy. The series consists of several manual aids which complete our wide range of patient lifts and accessories.*

*The Handy products also fulfil another need – the need to self participate according to one's ability. Especially for persons with disabilities, it is often very important to avoid becoming too passive and allow caregivers to do the entire job when being transferred in or outside the bed.*

*Several Handy products allow a high level of patient participation in the transfer and thereby support our philosophy not to use aids which unnecessarily make the patient passive.*

*First out in the series are the products presented in this brochure, HandySheet and HandyTube.*

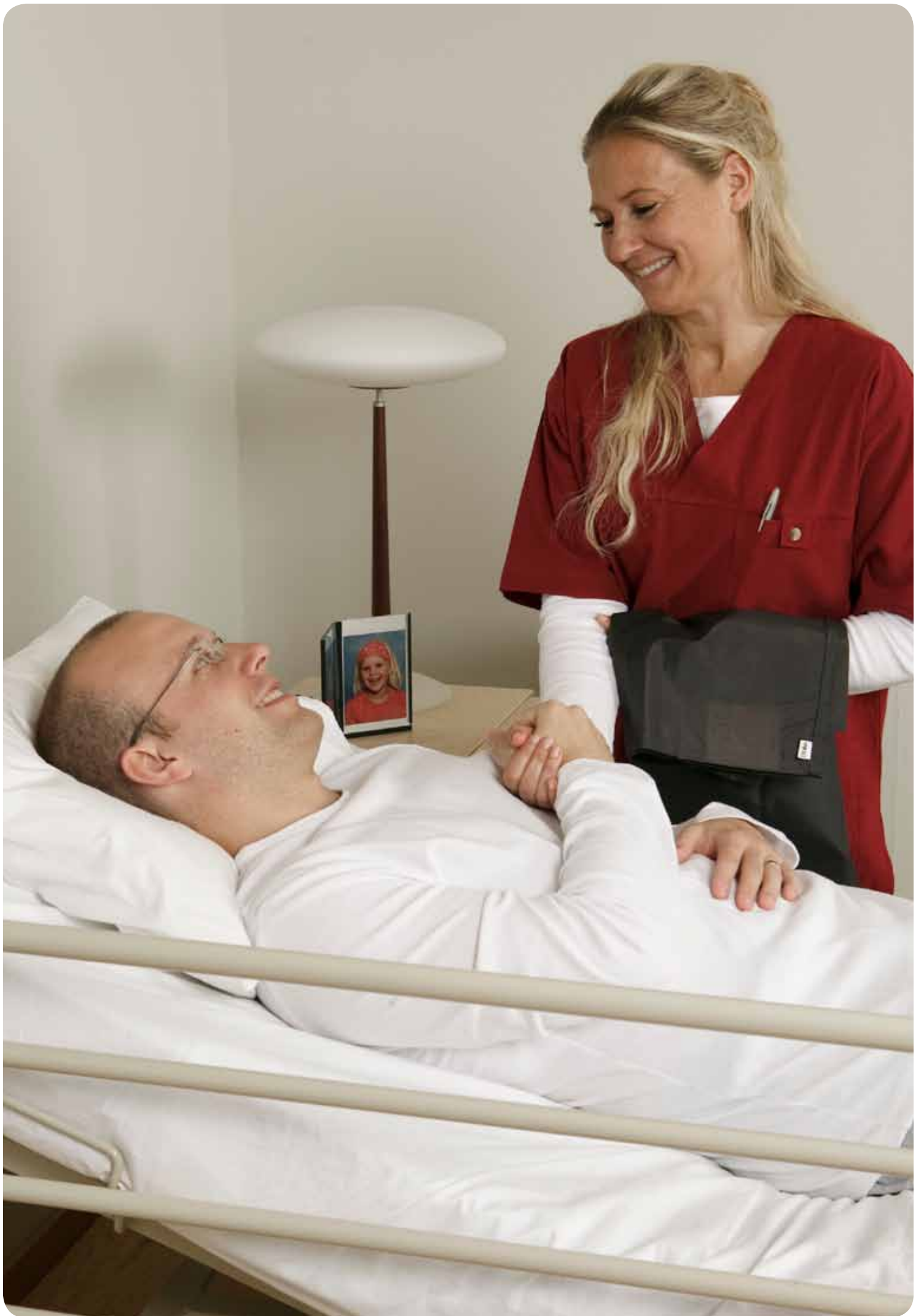



## HandySheet

HandySheet is an excellent aid to facilitate the application of slings, especially on bariatric patients, pain-sensitive patients and patients with contractures.

HandySheet is available in two sizes, the smaller model Short being mainly used to facilitate the application of the sling in a wheelchair or under the patient's legs. This model has small pockets which facilitate application in, for example, a narrow wheelchair.

The larger model Long is also used to apply a sling under a patient lying down. The HandySheet Long has loops in the upper edge to facilitate the attachment of the sliding sheet to the bed, thereby handling the operation with just two caregivers.

The sliding sheets are always delivered in a set of two, usually used in pairs.

Remember that a transfer with HandySheet should be simple and smooth. If the operations seem heavy, we recommend using other lift and transfer aids from our wide range.



### Applying a Sling on a Patient in Bed



1. For a patient lying down, two HandySheets Long are used.

2. The sliding sheets can be attached to the head end of the bed, using the smart loops and a FixStrap (accessory).

3. The two sliding sheets interact with each other and make it easy to insert and remove the sling from under the patient.



### Applying a Sling on a Sitting Patient



For a patient in a sitting position, the HandySheet Short makes it easier to insert the sling behind the patient's back ...



...or get the sling's leg supports into place under the patient's thighs.

#### PRODUCT DESIGNATION

	PROD. NO.	SIZE	MEASUREMENTS
HandySheet Short, without handles	3715015	M	70 x 70 cm (28 x 28 in.)
HandySheet Short, without handles	3715017	XL	120 x 85 cm (45 x 34 in.)
HandySheet Long, with handles	3716015	M	90 x 185 cm (36 x 74 in.)
HandySheet Long, with handles	3716017	XL	140 x 215 cm (56 x 86 in.)



# HandyTube

HandyTube is available in three different models with different application areas, and is packed individually (contrary to the HandySheet).

The smallest model, Short, is primarily used to move a patient higher up in the bed.

The medium model, Long, is used to turn a patient in the bed, or to transfer the patient between two beds.

HandyTube Wide can be used to turn or move a patient in the bed. It can also be used for transfers between two beds.

Remember that a transfer with HandySheet should be simple and smooth. If the operations seem heavy, we recommend using other lift and transfer aids from our wide range.

## PRODUCT DESIGNATION

	PROD. NO.	SIZE	MEASUREMENTS
HandyTube, Short	3710015	M	70 x 90 cm (28 x 36 in.)
HandyTube, Long	3711015	M	90 x 185 cm (36 x 74 in.)
HandyTube, Wide (with handles)	3712015	M	140 x 140 cm (56 x 56 in.)

### Turning in Bed



With HandyTube, the patient can be turned in bed more easily.



The patient is turned gently and in a controlled manner, using the low friction of the HandyTube.

### Higher Up in Bed



The patient can move higher up in bed on his/her own. The caregiver helps by holding the patient's feet.



The HandyTube enables repositioning with less effort, and the patient can participate actively.

### Between Two Beds



HandyTube Wide can be used to quickly and smoothly move a patient between two beds.



The need to move patients to and from transport stretchers is great, especially at hospitals.



### Liko has the complete solution for safe lifts and transfers

Now you have read about products which are an important part in solving lift and transfer situations in your business. But, actually, it is not enough. Over the years, we have found that a functioning and safe solution stands on three legs.

First, you need safe and easy-to-use **PRODUCTS**. Then, you need **KNOWLEDGE** on how to safely handle your aids and perform safe lifts and transfers. And finally, you need **WORK METHODS** which ensure that staff really work in the correct way – a program for safe transfers.

Through our knowledge, we provide the competence for all three parts. With Liko, you ensure that the staff as well as the patient are provided with a safe environment for lifting and transferring.

[www.liko.com](http://www.liko.com)

Manufactured by:  
Liko AB  
Nedre vägen 100  
975 92 Luleå  
Sweden  
[info@liko.se](mailto:info@liko.se)

Local representative in UK  
Hill-Rom Limited  
Clinitron House, Ashby Park  
Ashby de la Zouch,  
Leicester LE65 1JG  
Tel: 01530 411000  
Fax: 01530 411555  
[info@liko.co.uk](mailto:info@liko.co.uk)

To find your local Liko representative,  
please visit [www.liko.com](http://www.liko.com) or call  
+46 (0)920 47 47 00

 **Liko**  
A  Company